



INSTRUCTIONS FOR USE (IFU)

Fuse Plating Systems

Indications For Use:

The intended use of the Fuse Plating System is to draw two or more aligned bone fragments together in order to facilitate healing. It is composed of the following bone plate categories:

I. Small Bone Systems:

The Fuse Small Bone Plating System is indicated for fixation of fractures, osteotomies, non-unions, replantations, and fusions of short bones and short bone fragments including, but are not limited to, the hand, wrist, foot and ankle. The Fuse Small Bone Plating System is not for Spinal Use.

II. Foot System:

The Fuse FPS Foot Plating System is indicated for fixation of small bones and small bone fragments in the foot (Phalanges and Metatarsals), and medium/ large bones and medium/large bone multi-fragments in the foot (Cuneiform, Cuboid, Navicular, Talus and Calcaneus) for stabilization of fractures, joint fusions, osteotomies, nonunions, malunions, reconstruction of small, medium and large bones, revision surgeries and replantations in an adult patient. The Fuse FPS Foot Plating System is not for Spinal Use.

Contraindications Include:

- Infection.
- Patient conditions including blood supply limitations, obesity and insufficient quantity or quality of bone.
- Patients with mental or neurologic conditions who are unwilling or incapable of following postoperative care instructions.
- Foreign body sensitivity. If material sensitivity is suspected, testing is required prior to implanting the device.

R_x

Federal law restricts this device to sale by or on the order of a physician.

Materials:

The Fuse Plating System plates and screws are manufactured from a Titanium alloy (ASTM F136). The instruments are made of surgical grade stainless steel (ISO 7153-1 and ASTM F899).

Adverse Effects:

In all surgical procedures, the potential for complications and adverse reactions exist. The risks and complications with these implants include:


- Fracture of the implant due to excessive loading
- Incomplete or inadequate healing
- Implant migration and / or loosening
- Pain, discomfort, or abnormal sensations due to the presence of an implant
- Nerve damage resulting from surgical trauma.
- Bone necrosis or bone resorption
- Delayed or nonunion of bone fragments
- Allergic reaction to the implant materials

Warnings & Precautions:

- Re-operation to remove or replace implants may be required at any time due to medical reasons or device failure. If corrective action is not taken, complications may occur.
- Implants must not be re-used or re-sterilized.
- Improper insertion of the device during implantation may result in implant loosening or migration.
- Loosening or migration and loss of fixation due to incorrect implantation, delayed union, nonunion and incomplete healing may occur.
- Bending or fracture due to applied excessive stresses and load bearing.
- Failure to follow postoperative care instructions may result in procedure complications or failure.
- Electrolytic action and corrosion due to implanting with other metallic devices of different chemical composition may occur.

MRI Safety Information:

MRI Safety Information



A patient with the Fuse Plating System may be safely scanned under the following conditions. Failure to follow these conditions may result in injury to the patient.

Device Name/Identification	Fuse Plating System
Nominal value(s) of Static Magnetic Field (T)	1.5 T or 3.0 T
Maximum Spatial Field Gradient [T/m and gauss/cm]	20 T/m (2000gauss/cm)
RF Excitation	Circularly Polarized (CP)
RF Transmit Coil Type	Whole body transmit coil, Head RF transmit-received coil
Maximum Whole Body SAR [W/kg]	1.0 W/kg or 2.0W/kg (Normal Operating Mode)
Limits on Scan Duration 1.0 W/kg SAR	1.0 W/kg whole body average SAR for 60 minutes of continuous RF (a sequence or back-to-back series/scans without breaks)
MR Image Artifact	The presence of this implant may produce an image artifact of 22 mm.

Instructions For Use:

1. Using standard dissection techniques, expose the surgical site.
2. Perform the intended osteotomy or identify the fracture location.
3. After reduction of the fracture, choose the proper plate based on the size and type of indication.
4. Place the plate on the fracture/osteotomy site, fix with k-wires. If forming/bending the plate to fit the anatomy – use the bending irons for preparation of the proper contour. **DO NOT REPEATEDLY BEND THE PLATE** – as this will cause a weakened fatigue life of the plate.
5. Utilize the drill guide with proper drill according to screws size for angulation into the most secure bone structure. Drill hole for screw. Repeat hole preparation as necessary for proper fixation of the plate.
6. Utilize the depth gauge for proper length of screw in bone anatomy for firm fixation in the opposite bone cortex.
7. Insert desired size screw matching to plate size and bone anatomy. Repeat process on remaining screw(s) with angulation holes using either locking or non-locking screws.
8. Remove k-wires. Check plate/screw tightness on bone anatomy fracture/osteotomy site.

9. Using fluoroscopy, confirm the proper plate and screw placement on the bone anatomy. Correct as warranted & re-check.
10. Clean the surrounding area with a pulse lavage.
11. Use the surgeon's preferred method for closing the surgical site.

Postoperative Management:

The patient is allowed to ambulate with weightbearing to tolerance on the operated fracture site within limits imposed by postoperative discomfort. The progression to normal use of the digit or limb is limited only by the persistence of postoperative swelling and discomfort.

Care and Handling:

Certain components are provided non-sterile and should be stored in the original packaging until cleaned and sterilized. Prior to use, they must be sterilized according to the standard hospital procedure. Refer to the STERILIZATION section for recommended parameters.

Limitations on Processing:

Repeated processing has minimal effect on these implant and instruments. End of life is normally determined by wear and damage due to use.

Point of Use:

Before being used for the first time and each use thereafter, the instructions outlined below should be followed to ensure safe handling of biologically contaminated instruments.

Containment and Transportation:

It is recommended that instruments are reprocessed as soon as reasonably practical following use.

Preparation for Cleaning

Where instruments interface with other devices, disassemble prior to cleaning. Remove excess soil with a clean, disposable, absorbent Kimwipe or equivalent.

Disassembling of Depth Indicator:

1. Unthread proximal cap counter- clockwise until cap disengages from outer cannula of instrument.
2. Remove cap and inner cannulas from outer cannula.
3. Proceed to cleaning steps below.

Cleaning (Automated):

Equipment: Automated washer, soft bristle brush, enzymatic detergent, and neutral pH detergent.

- Preclean the instruments by placing them under running water and scrubbing with a soft bristle brush to

remove major debris. Rinse and scrub each instrument for at least one minute.

- After precleaning, place in the automated washer, making sure the samples do not touch each other - load instruments in such a way that the parts can drain.
- Use a standard instruments cycle with the following parameters (at a minimum):

Enzyme Wash	Hot 40-65 °C / 104-149 °F for 3 minutes
Neutral pH Wash	60 °C / 140 °F for 3 minutes
Rinse	Ambient temperature for 1.5 minutes
Thermal Rinse	90 °C / 194 °F for 1 minute
Dry	82 °C / 180 °F for 6 minutes

- Determine if the instruments are dry. If they are not dry, dry with a soft, clean, lint free cloth.
- After drying, check instruments for complete removal of any debris. If necessary, repeat cycle or use manual cleaning.
- Final Rinse shall be performed, using reverse osmosis or distilled water at room temperature for 5 minutes.
- The rinse bath should be changed after each cleaning process.

Cleaning (Manual):

Warning: Movable components and blind holes require particular attention during cleaning.

Preparation of Cleaning Agents (Recommended):

- Add 60 mL of Endozime® AW Plus to 3.8 L of water, (1:64 dilution).

Manual Cleaning Instructions:

- Preclean the instruments by placing them under running water and scrubbing with a soft bristle brush to remove major debris. Rinse and scrub each instrument for at least one minute.
- Bathe the instruments in the enzymatic solution for 5 minutes; where appropriate, the instrument shall be rotated and briskly moved in bath to promote flushing. Where appropriate, a large syringe or pulsating water jet may be used to thoroughly flush all channels and lumens with the solution.

- Scrub the instruments with a soft bristle brush while submerged in the detergent.
- Rinse the devices using reverse osmosis or distilled water at room temperature for 5 minutes.
- The rinse bath should be changed after each cleaning process.
- Pat dry with a soft, clean, lint free cloth.
- After drying, check instruments for complete removal of any debris. If necessary, repeat manual cleaning.

After Cleaning:

Where instruments have been disassembled prior to cleaning reassemble prior to use.

Inspection and Function Testing All instruments:

Visually inspect for damage and wear. Where instruments interface with other devices, inspect to ensure that the interface is not damaged. Check for misalignment, burrs, bent or fractured tips. Mechanically test the working parts to verify that each instrument functions correctly. Remove stained, discolored, or damaged instruments.

Packaging:

Instruments may be loaded into the specified instrument trays, or general-purpose trays. Wrap the trays using an appropriate method with no more than two layers of sterilization wrap that are FDA cleared for pre-vacuum steam sterilization.

Sterilization:

For components provided Sterile, the sterilization method is noted on label. Sterile implant components are supplied sterile to a Sterility Assurance Level (SAL) of 10⁻⁶. Sterile packaged components are supplied in a protective sterile barrier packaging. Inspect package for punctures or other damage prior to surgery. If the sterile barrier has been broken, return component to Fuse Medical.

WARNING: Please note that a single-use device (SUD) which comes in contact with human blood or tissue should not be re-used and should be returned to the manufacturer or properly disposed.

If not specifically labeled STERILE, or if labeled NON-STERILE, components are supplied non-sterile and must be cleaned and sterilized prior to surgery.

Warning: It is not recommended that the instruments be sterilized by Flash, EtO or Chemical sterilization. When sterilizing multiple instruments in one autoclave

cycle, ensure that the sterilizer's maximum load is not exceeded.

To achieve a sterility assurance level of SAL 10⁻⁶, the following parameters are recommended:

Method	Time	Temperature	Dry Time
Pre-vacuum	4 minutes	270° F (132° C)	20 minutes
	3 minutes	275° F (135° C)	
Gravity	15 minutes	270° F (132° C)	20 minutes

It is recommended to follow ANSI/AAMI ST79, Comprehensive guide to steam sterilization and sterility assurance in health care facilities, which includes: physical monitoring of the cycle, inclusion of a chemical indicator internal and external to the package, and monitoring of every load with a Biological Indicator and/or Class 5 Integrating Indicator.



Storage:





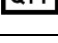


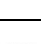
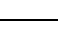




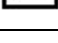
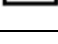
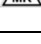
Fuse Plating System instruments must be completely dry before storing and must be handled with care to prevent damage. Store in designated trays and in areas which provide protection from dust, insects, chemical vapors and extreme changes in temperature and humidity.

Retrieval and Analysis of Removed Implants:

The most important part of surgical implant retrieval is preventing damage that would render scientific examination useless. Special care should be given to protect the implant during handling and shipping. Follow internal hospital procedures for the retrieval and analysis of implants removed during surgery. When handling removed implants, use precautions to prevent the spread of bloodborne pathogens. Please contact Fuse Medical customer service prior to the return of removed implants.

Symbol Glossary

Symbol	Title/ Standard	Meaning
	21 CFR 801.109b Prescription Only	Indicates that a practitioner licensed by the law of the state in which the practitioner practices to use or order the use of the device
	ISO 15223-15.1.6 Catalogue Number	Indicates the manufacturers catalogue number so that the medical device can be identified.

	ISO 15223-15.1.5 Batch Code	Indicates the manufacturer's batch code so that the batch or lot can be identified
	Material	Indicates the material of the device
	ISO 15223-1 5.1.11 Country of manufacture	To identify the country of manufacture of products
	ISO 15223-15.1.4 Use-by date	Indicates the date after which the medical device is not to be used
	Quantity	Indicates the quantity of devices.
	ISO 15223-15.2.4 Sterilized using irradiation	Indicates a medical device has been sterilized using irradiation.
	ISO 15223-15.4.2 Do not re-use	Indicates a medical device that is intended for one single use only.
	ISO 15223-15.2.8 Do not use if package is damaged	Indicates a product should not be used if the package is damaged or opened and the user should consult the instructions for use for additional information.
	ISO 15223-1 5.4.3 Consult instructions for use	Indicates the need for the user to consult the instructions for use.
	ISO 15223-1 5.4.3 Caution	Indicates the need for the user to consult the instructions for use for important cautionary information such as warnings and precautions that cannot, for a variety of reasons, be presented on the medical device itself.
	ISO 15223-1 5.2.7 Non-sterile	Indicates a medical device that has not been subjected to a sterilization process
	ISO 15223-1 5.1.9 Distributor	Indicates the entity distributing the medical device into the locale.
	ISO 15223-1 5.1.1 Manufacturer	Indicates the medical device manufacturer.
	ISO 15223-1 Medical Device Symbol	Indicates that the item is a medical device.
	ISO 15223-1 5.7.10 Unique Device Identifier	Indicates a carrier that contains unique device identifier information.
	ASTM F2503	Indicates that the device is MR Conditional and can be used in the MRI environment provided certain strict conditions are followed.

Additional Information:

Comments, questions, or concerns regarding the use of this device should be directed to the following:

Website
Email
Phone
Postal

www.FuseMedical.Com
Info@FuseMedical.Com
469-862-3030
Attn: Customer Service
CPM Medical Consultants, LLC
4343 Sigma Road, Suite 500
Farmers Branch, TX 75244
www.FuseMedical.Com

Website



Distributed By:
CPM Medical Consultants, LLC
4343 Sigma Rd., Suite 500
Farmers Branch, TX 75244



Manufactured By
Tyber Medical
83 S Commerce Way, Suite 310
Bethlehem, PA 18017