



## DERMAL MATRIX

**FuseChoice™ Derm and FuseChoice™ Derm X** are non-fenestrated, de-epithelialized, decellularized, & disinfected collagen matrices indicated for supplemental support, protection, reinforcement, or covering of a tendon or soft tissue.

### COMPOSITION

- Acellular dermal matrix allograft with two thickness options to provide optimal integumentary augmentation support.
- Applied as reinforcement to help provide protection and framework while retaining intact vascular channels.<sup>1</sup>
- Proprietary process preserves the inherent properties of dermal matrices, maintaining inherent levels of key extracellular matrix molecules, including collagen and elastin.<sup>1</sup>



### VERSATILITY

- Excellent tensile strength to maintain framework.<sup>1</sup>
- Biocompatibility study demonstrated attachment and proliferation of:
  - Fibroblast cells
  - MIAMI cells
  - Primary chondrocytes
- Variety of sizes and thicknesses available

### SAFETY

- Requires no up-front preparation.
- Unique processing technique preserves the collagen and elastic tissue fibers while maintaining the open channels through which cells can migrate, proliferate, and form new blood vessels.
- Ambient temperature storage with a 5-year shelf life.
- E-Beam sterilization provides sterility assurance level (SAL) of  $10^{-6}$ .

## DERMAL MATRIX THIN 0.4 - 0.8mm

- Ambient Temperature (2°C to 30°C)
- 5-year shelf life
- Non-fenestrated
- Sterile packaged
- Affix using a minimum of 2-0 suture

### FuseChoice™ Derm


PRODUCT NUMBER	SIZE
FCD-5440	4x4 cm
FCD-5480	4x8 cm

## DERMAL MATRIX THICK 1.0 - 2.0mm

- Ambient Temperature (2°C to 30°C)
- 5-year shelf life
- Non-fenestrated
- Sterile packaged
- Affix using a minimum of 2-0 suture

### FuseChoice™ DermX

PRODUCT NUMBER	SIZE
FCDX-5440	4x4 cm
FCDX-5480	4x8 cm

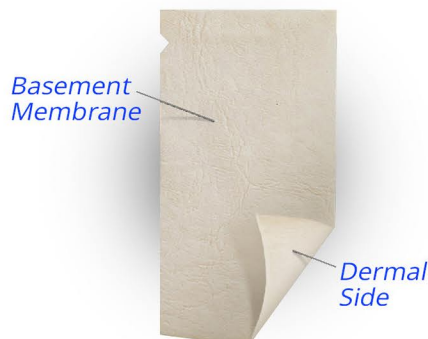

To orient the graft with the basement membrane facing up, the notch must be in the upper left hand corner of the graft (see figure to the left).

**fuse**

1. Data on file at VIVEX Biologics, Inc.

Fuse Medical, Inc. will use reasonable efforts to provide accurate and complete information herein, but this information should not be construed as providing clinical advice, dictating reimbursement policy or as a substitute for the judgment of a healthcare provider. It is the healthcare provider's responsibility to determine the appropriate treatment, codes, charges for services and use of modifiers for services rendered and to submit coverage or reimbursement-related documentation.

**Fuse Medical, Inc.**

1565 N. Central Expressway, Suite 220  
Richardson, TX 75080  
469.862.3030

info@fusemedical.com  
www.fusemedical.com

DC-0116