


# **SILKTOE**



Metatarsophalangeal Joint Arthroplasty System

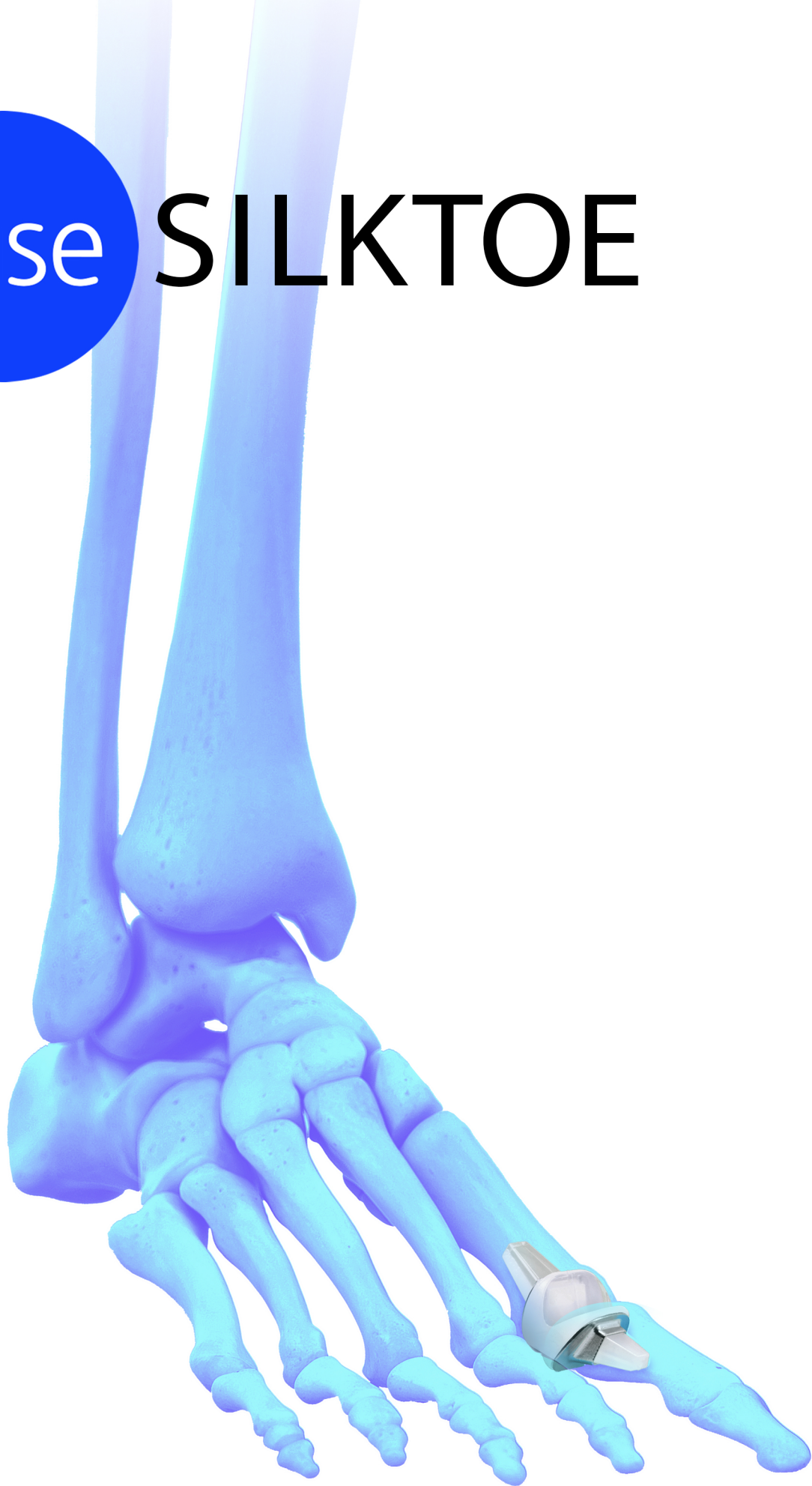


## System Highlights

Silktoe® metatarsophalangeal joint arthroplasty system offers:

- Right and Left specific options
- Elastomeric medical-grade silicone
- 4-degrees of dorsiflexion built-in to restore anatomical motion
- Distal stem - trapezoidal shape contours to the intramedullary canal of the phalanx
- Proximal Stem - rectangular shape contours to the intramedullary metatarsal canal
- Sterile packaged
- 20, 30, 35, 40, & 50mm length options

### Ring Specifications

- Titanium
- Acts as a support stabilizer for the silicone joint spacer
- Sterile packaged

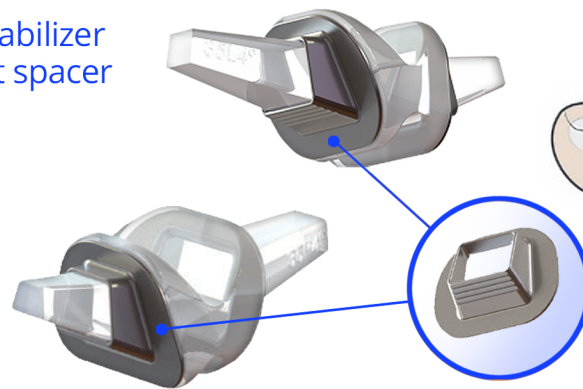

**The intended use of the Silktoe® device is to represent a permanent implant, used in all severe forms of arthritis, such as:**

- Hallux rigidus or hallux limitus
- Painful rheumatoid arthritis
- Hallux abducto valgus associated with arthritis
- Unstable or painful joint from previous surgery

Fuse Medical, Inc.

4343 Sigma Rd., Suite 500  
Farmers Branch, TX 75244  
469.862.3030

info@fusemedical.com  
www.fusemedical.com

SILTB-001 REV B DC-0467